

Safety Precautions for Product Use


The Lift2Glow LIFT Cartridge is a high-power, energy-based device.
As its effects are significant, there is a risk of skin damage if the instructions are not strictly followed.

Why is the LIFT Cartridge handled with care? Lift2Glow is designed to deliver energy deep into the skin to improve elasticity. Therefore, if the device is not in full contact with the skin or if used at a high level from the start, localized heat stimulation may concentrate and lead to skin damage.

4 Essential Safety Rules



- 1. Full Skin Contact:** Only press the SHOT button when the cartridge head is in complete contact with the skin.
- 2. Start at "Level 1":** When using for the first time, you must start at Level 1. As skin sensitivity varies by individual, increase the intensity step-by-step.
- 3. Apply Gel:** You must use a **Device Compatible (Aloe Type)** gel during treatment.
- 4. Keep Moving:** Do not stay on one spot for too long or treat the same area repeatedly in a single session.

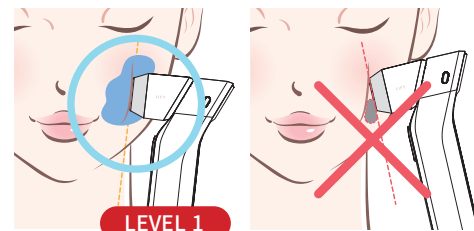

10-Second Safety Check Before Use!


The key to great results is safe usage.

- ✔ Have you applied a sufficiently thick layer of **Device Compatible (Aloe Type)** gel?
- ✔ Is the LIFT cartridge head in full contact with your skin?
- ✔ Is this your first time using it today? **Start at Level 1!**
- ✔ **Avoid these areas:** Between eyebrows, eye area, eyeballs, philtrum, lips, around the mouth, nose, jawbone, Adam's apple line, and thyroid area.

SHURINK
— HOME —

LIFT **2** GLOW Quick Start Guide

How to Adjust Energy Level and Mode

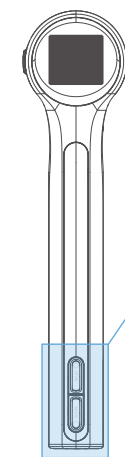


Voice Setting

: Press and hold for 2 seconds to cycle through the following voice modes in order.
(Korean→English→Mute)


Energy Level

: Press once.

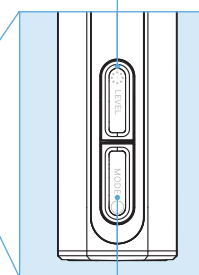

Power/Mode Control

- Press and hold: Turn ON/OFF
- Press once: Switch Mode

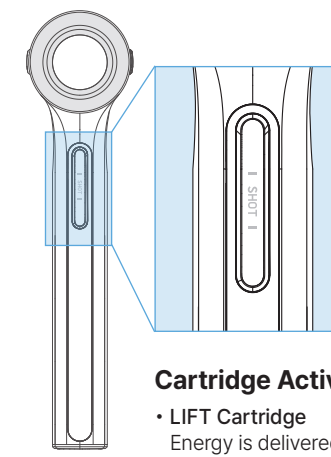


Cartridge Activation

- **LIFT Cartridge**
Energy is delivered while the button is being pressed and held.
- **GLOW Cartridge**
Press once to deliver a single energy shot.

Display Guide

LIFT Cartridge

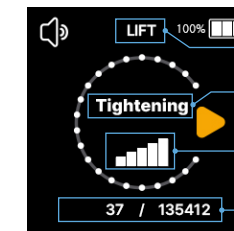

- Cartridge Type
- Care Mode
- Energy Level
- Current Shot Count / Total Shot Count

GLOW Cartridge



- Cartridge Type
- Energy Level
- Care Duration (Minutes : Seconds)

LIFT Cartridge Care Mode User Guide

Contour



Gently place the device on your skin, then roll slowly and smoothly from bottom to top along the jawline, cheeks, and cheekbones for 5 minutes of care.


Tightening



Gently place the device on your skin, then lift slowly along the jawline, cheeks, and cheekbones for 5 minutes of care.

*** For optimal results, move the device gradually to the next area each time you hear the sound indicating a single shot.**

Eye/Neck

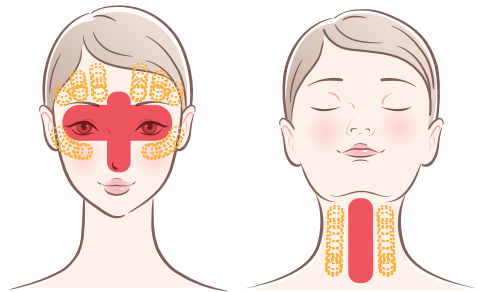

Gently place the device on your skin, then slowly lift and care for the areas around the eyes and forehead for 3 minutes each, and the neck area for 2 minutes.

(Warning! Do not use on areas marked in red)

Smile Line



Gently place the device on your skin, then slowly roll it around the mouth area and gently care for 2 minutes.

GLOW Cartridge Care Mode User Guide

Total Care Duration: 6 minutes (3 minutes on each side: Left and Right)




Jawline

- Gently move from the center of the chin toward the ears, gently rolling upward from bottom to top.
- To enhance circulation, roll again over the lymph nodes.

Cheekbones

- Gently roll upward along the cheekbones, starting from the corners of the mouth.

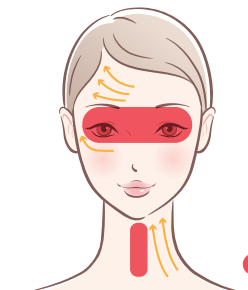

Warning
: Avoid These Areas

Eyes

- Gently roll under the eyes from the inner to the outer corner.
- Avoid using on the eye area.

Forehead


- Gently roll upward from the center of the forehead toward the temples.

Neck

- Gently roll upward from the collarbone to the jawline.
- Avoid using on the throat (Adam's apple) and thyroid area.

LIFT Cartridge Warning: Avoid These Areas

- Do not use on the area between the eyebrows, eye area, lips, nasal bone, jawline, throat line, and thyroid.
- The manufacturer is not responsible for any side effects resulting from use on the restricted areas.

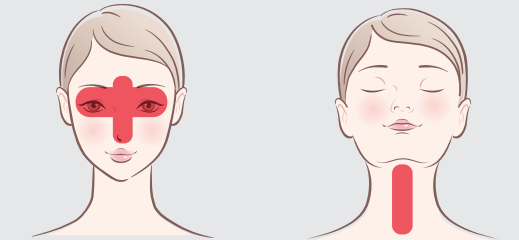

Safety Precautions

LIFT Cartridge

- To avoid the risk of burns, continuously move the device and shift to a new area whenever you hear the single-shot notification sound.
- Please evenly apply the gel cream. Insufficient application of gel cream may reduce the effectiveness of energy delivery.

GLOW Cartridge

- Before use, apply a generous amount of gel and ensure the device maintains full contact with the skin for continuous operation.
- Insufficient gel or poor contact with the skin may cause ultrasound waves to reflect, leading to excessive heat and increasing the risk of burns.
- When removing the device from the skin during use, always press the shot button to temporarily pause operation before lifting it.